
Gas Transmission Infrared Fiber Optic Probe for Gases



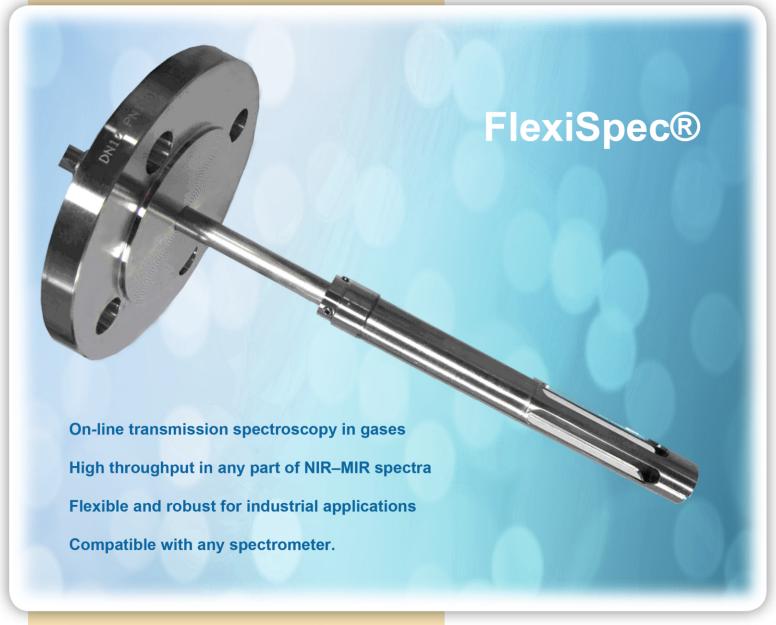

FlexiSpec® product line includes the latest generation of Transmission infrared fiber optic probes for gases to be used with any spectrometer or photometer. The design of Transmission Fiber Probe for gases is based on a bifurcated fiber bundle.


Due to low Mid IR-attenuation in gases the collimated beam design of double pass (or multiple pass) gas cell is needed to increase optical path length to 10 – 40cm. This design is easily realized by the mean of the collimating objective and reflecting mirror cell.

Applications:

- Reaction investigation and monitoring in real time
- Analytical Characterization
- Exhaust gases monitoring
- Solvent vaporization monitoring
- Associated petroleum gas monitoring

broad spectra fiber solutions


www.artphotonics.com

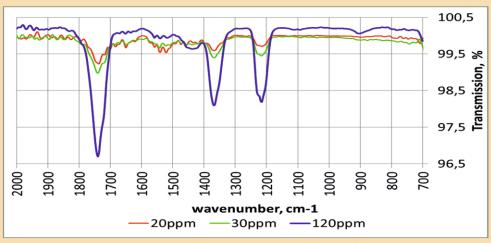

Acetone vapor spectra in 100mm double-pass cell

Specification of Gas Transmission Fiber Optic Probes FlexiSpec®				
Probe type	Transmission in double-pass gas cell			
Transmission range	0,5 - 2,2 μm	1,6 – 5,5 μm	3 – 18µm	
Fiber type	Silica	Chalcogenide IR	Polycrystalline IR	
Temperature range	-50°C + 200°C	-50°C + 90°C	-50°C + 140°C	
Sensitivity	>1ppm depending on optical path, spectrometer performance, gases in test			

Common Parameters of Gas Transmission Fiber Optic Probes FlexiSpec®		
Total Length	1,5 m (opt: 1 – 30m) *	
Transmission Cell Length	100mm (opt: 5 – 200mm) *	
Transmission cell Diameter	25 mm	
Shaft Material	Stainless Steel, Hastelloy C22	
Protective Tube Material	Liquid Tight SS-Conduit, KOPEX-Tube	
Input / Output Connectors	Long SMA *	
*Customized dimensions are available on request		

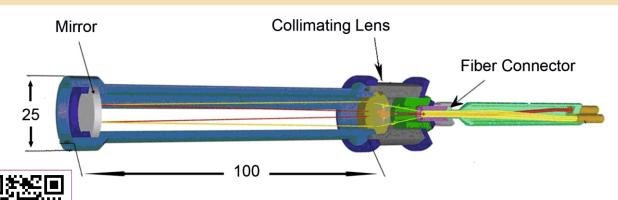

Phone + 49 (0) 30 6789 4153

sales@artphotonics.com

www.artphotonics.com



via smartphone